


# DART SERVICE INSTRUCTION

TO AMEND INSTALLATION INSTRUCTIONS IIN-D350-616 REV. C


REF TCCA STC SH96-10

REF FAA STC SR00405NY

IF REPLACING THE D2493 PATIENT STOP, IT MAY BE NECESSARY TO TRANSFER DRILL THE PIP PIN ATTACHMENT HOLE FROM THE D2370 LITTER ASSEMBLY TO THE D2493 PATIENT STOP. DRILL Ø0.203" (5.2mm) TO MATCH HOLE IN LITTER ASSEMBLY. TOUCH-UP PAINT AS REQUIRED TO MATCH ORIGINAL FINISH.

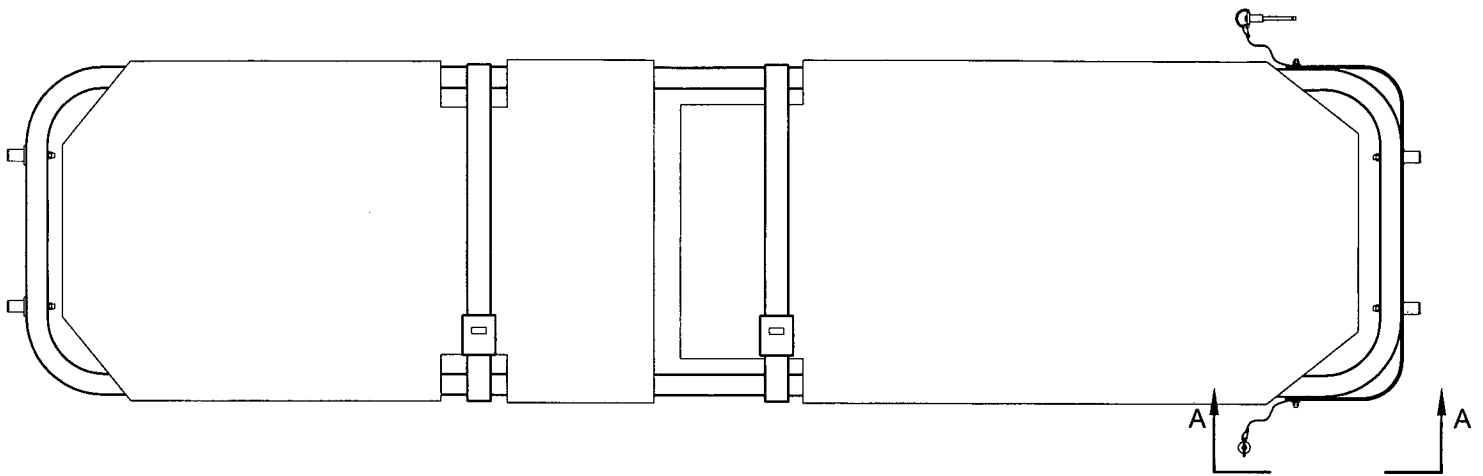

**FIGURE 1: LITTER ASSEMBLY**  
(VIEW LOOKING DOWN)

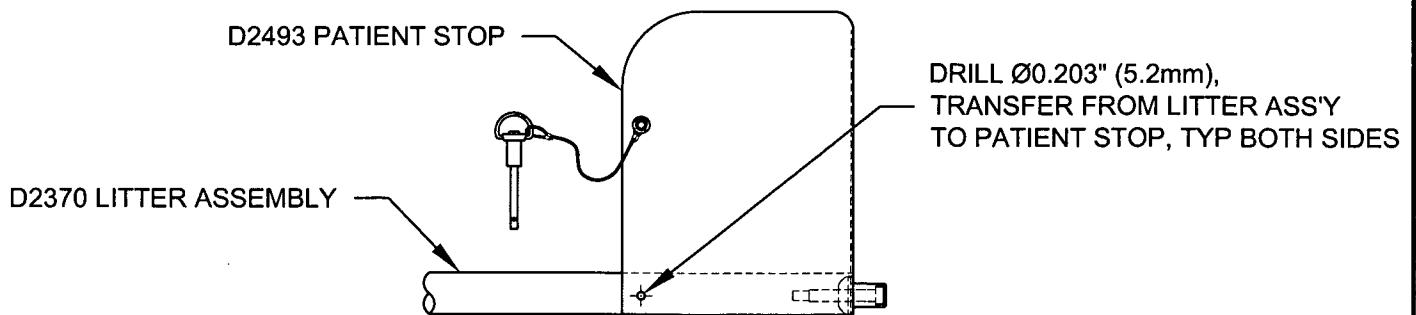

**VIEW A-A**  
(SCALE 2X)

CANADA  
DEPARTMENT OF TRANSPORT  
AIRCRAFT CERTIFICATION  
BRANCH  
DAO # 01-O-01

APPROVED  
BY: *[Signature]*  
D. SHEPHERD (DE # 02)

DATE: 10.08.12  
CERT. NO.: SH96-10  
ISSUE NO.: 2

A	NEW ISSUE (NCR 10-062)		CP	10.08.12
REV.	DESCRIPTION		BY	DATE
DESIGN	99	<b>DART AEROSPACE LTD</b> HAWKESBURY, ONTARIO, CANADA		
DRAWN	99			
CHECKED	N/A	DRAWING NO.	REV. A	
MFG. APPR.		DSI 9530	SHEET 1 OF 1	
APPROVED		TITLE	SCALE	
DE APPR.		PATIENT STOP REPLACEMENT	NTS	
DATE	10.08.12		COPYRIGHT © 2010 BY DART AEROSPACE LTD	
THIS DOCUMENT IS PRIVATE AND CONFIDENTIAL AND IS SUPPLIED ON THE EXPRESS CONDITION THAT IT IS NOT TO BE USED FOR ANY PURPOSE OR COPIED OR COMMUNICATED TO ANY OTHER PERSON WITHOUT WRITTEN PERMISSION FROM DART AEROSPACE LTD.				